


# Certificate of Analysis

V<sub>1</sub>-1

ISO GUIDE 34  
ANAB Cert# AR-1470

ISO/IEC 17025  
ANAB Cert# AT-1467

## URACIL CERTIFIED REFERENCE MATERIAL

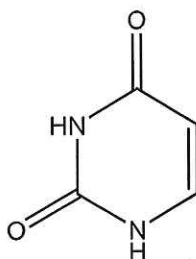

**CERTIFIED PURITY: 99.9%,**  $U_{\text{crm}} = \pm 0.05\% \text{ } k = 2$   
(Mass Balance/as is basis)

**NOMINAL PACKAGE SIZE: 100MG**

**CATALOG #:** PHR1581

**LOT #:** LRAB3029

**CERTIFICATE VERSION:** LRAB3029.2

**ISSUE DATE:** 30 April, 2017

*Note: Certificates may be updated due to Pharmacopeial Lot changes or the availability of new data.*

*Check our website at: [www.sigma-aldrich.com](http://www.sigma-aldrich.com) for the most current version.*

**CRM EXPIRATION:** 31 December 2021 (Proper Storage and Handling Required).

**RECEIPT DATE:** \_\_\_\_\_

Note: this space is provided for convenience only and its use is not required.

**STORAGE:** Store in a Refrigerator/Protect from Light, keep container tightly closed.  
Attachment of a 20 mm aluminum crimp seal recommended for unused portions.

**CHEMICAL FORMULA:** C<sub>4</sub>H<sub>4</sub>N<sub>2</sub>O<sub>2</sub>

**MW:** 112.1

**PHYSICAL DESCRIPTION:** White powder in amber vial **CAS #:** 66-22-8

**HAZARDS:** Read Safety Data Sheet before using. All chemical reference materials should be considered potentially hazardous and should be used only by qualified laboratory personnel.

**SIGMA-ALDRICH®**



Page 1 of 6

**INSTRUCTIONS FOR USE:** Do not dry, use on the as is basis. The internal pressure of the container may be slightly different from the atmospheric pressure at the user's location. Open slowly and carefully to avoid dispersion of the material. This material is intended for R&D use only. Not for drug, household or other uses.

## **TRACEABILITY ASSAY**

Comparative assay demonstrates direct traceability to Pharmacopeial Standards

**METHOD:** HPLC (ref.: Adapted from Fluorouracil, Current Compendial Monographs)

### ***ASSAY vs. USP REFERENCE STANDARD (as is basis)***

**ASSAY VALUE**

**99.1%**

**vs. USP LOT**

**F1K302**

Labeled Content = 1.00mg/mg

Column: Ascentis Express C18, 4.6 x 100mm, 5µm

Mobile Phase: 6.805g/L KH<sub>2</sub>PO<sub>4</sub> in Water (pH: 5.7)


Flow Rate: 1mL/min

Column Temperature: 30°C

Injection: 10µL

Detector: 266nm

Representative Chromatogram from Lot: LRAB3029 Analysis

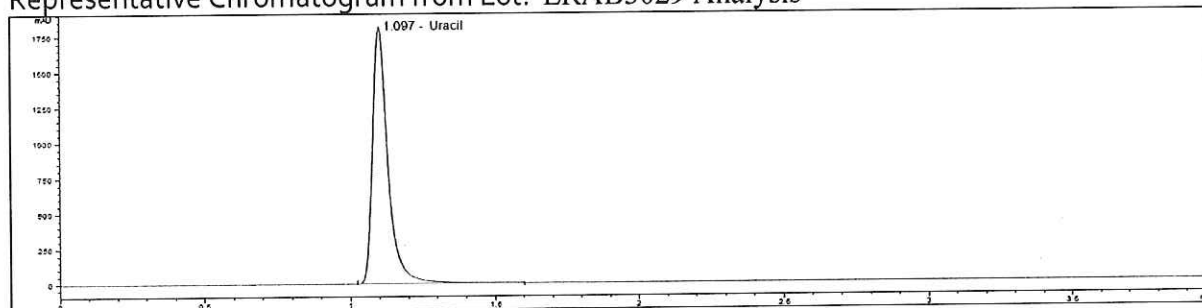

### ***ASSAY vs. EP CRS (as is basis)***

**ASSAY VALUE**

**100.4%**

**vs. EP BATCH**

**1.1**

Labeled Content = None

Assigned Content = 100.7%\*


\*The assigned content of the EP CRS was determined by assay against the USP Reference Standard

**SIGMA-ALDRICH®**



Column: Ascentis Express C18, 4.6 x 100mm, 2.7µm  
Mobile Phase: 6.8g/L KH<sub>2</sub>PO<sub>4</sub> in Water (pH: 5.7)  
Flow Rate: 1mL/min  
Column Temperature: 30°C  
Injection: 10µL  
Detector: 266nm

Representative Chromatogram from Lot: LRAB3029 Analysis

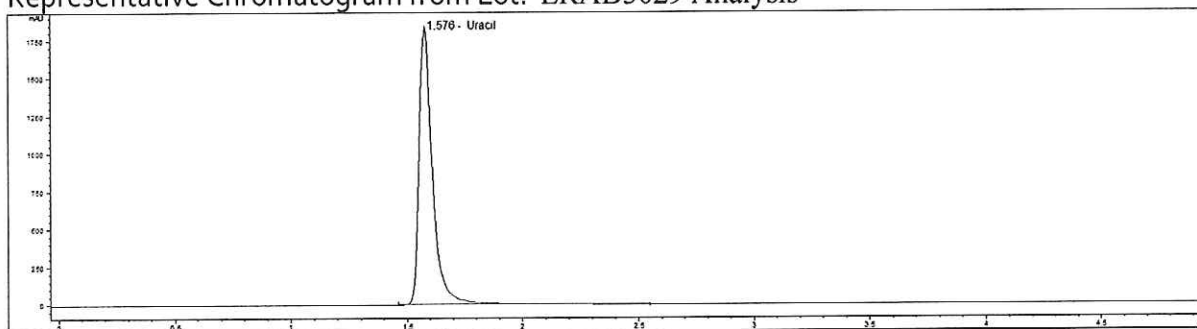

## PURITY DETERMINATION BY MASS BALANCE


### CHROMATOGRAPHIC IMPURITY ANALYSIS METHOD: HPLC

See Assay

Impurities Detected:

Impurity 1: 0.004%

Representative Chromatogram from Lot: LRAB3029 Impurity Analysis

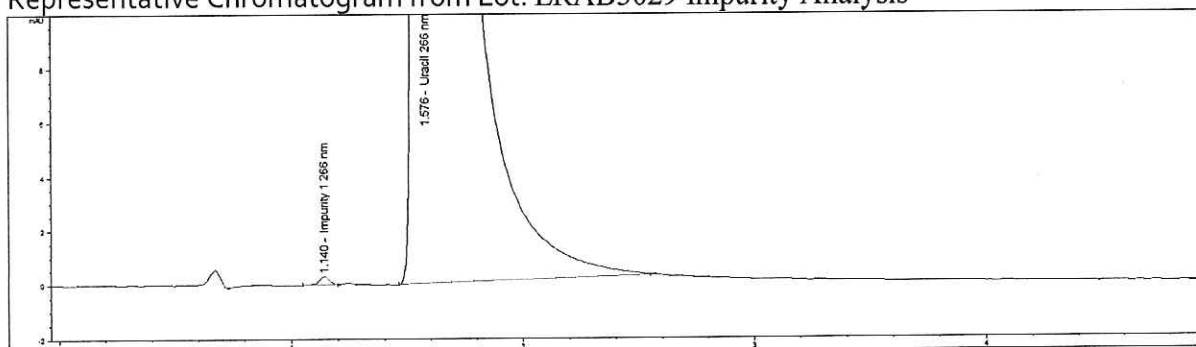

**SIGMA-ALDRICH**





## RESIDUAL SOLVENTS

Method: GC-MS Headspace (ref.: Adapted from Residual Solvents USP <467>)

Column: DB-1301

Carrier gas: He

Flow: 1.2mL/min

Split Ratio: 1:5

Injection/Temperature: 1mL/250°C

Temperature Program: 40°C for 20min, 10°C/min to 240°C, hold 20min

Solvents Detected: None

## LOSS ON DRYING/VOLATILES

Method: Oven at 105°C

Mean of three measurements, Loss = **0.05%**

## RESIDUE ANALYSIS

Method: Sulfated Ash

Sample Size: ~500mg


Mean of three measurements, Residue = **0.003%**

## CERTIFIED PURITY BY MASS BALANCE [100% - Impurities (normalized)]

**99.9%**  $U_{\text{crm}} = \pm 0.05\%$ ,  $k = 2$   
(as is basis)


## IDENTIFICATION TESTS

**INFRARED SPECTROPHOTOMETRY** (Comparative identification analysis demonstrates direct traceability to Pharmacopeial standards)

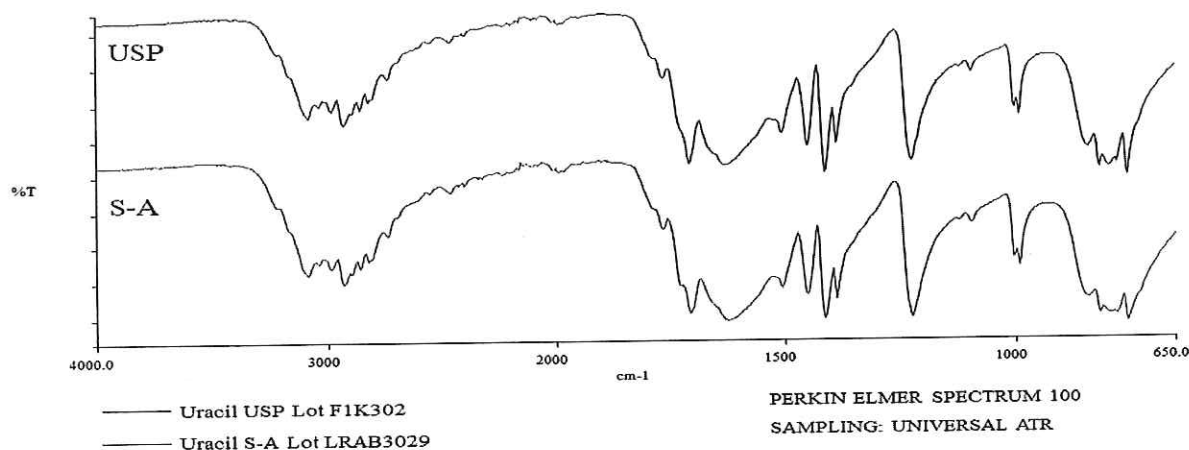

**SIGMA-ALDRICH®**



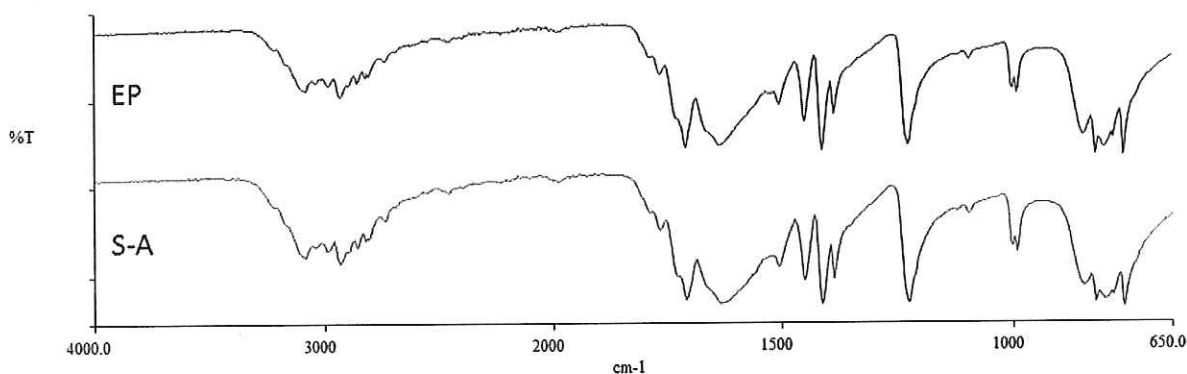


— Uracil EP Batch 1.1  
— Uracil Lot LRAB3029

PERKIN ELMER SPECTRUM 100

SAMPLING: UNIVERSAL ATR

**$^1\text{H}$  NMR** (Data provided by an external laboratory; not in scope of accreditation)

LRAB3029 Uracil

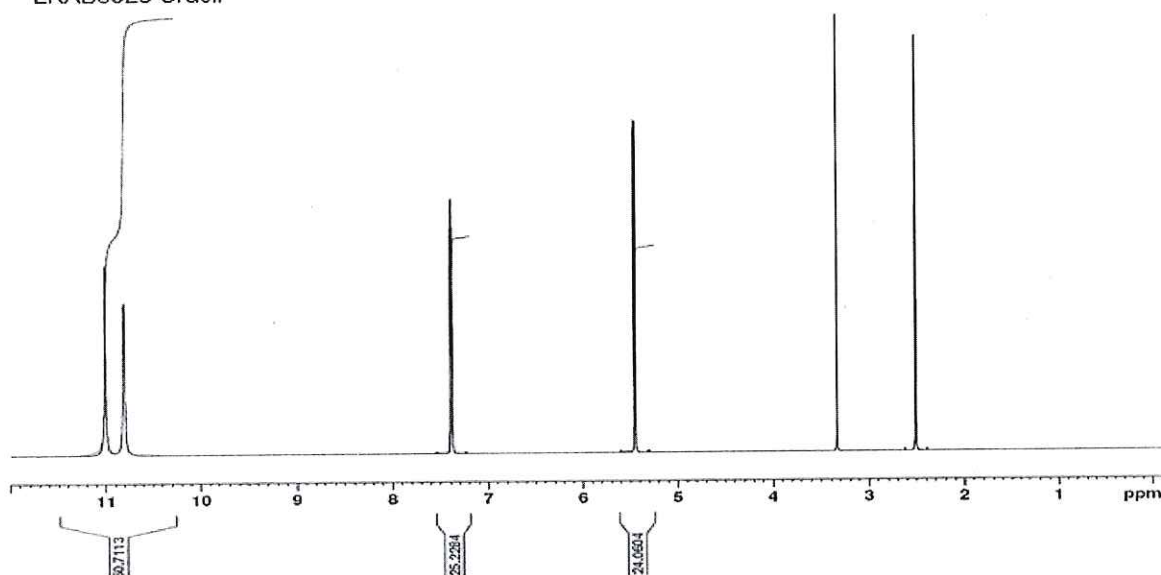

Consistent with structure

**ELEMENTAL ANALYSIS** (Data provided by an external laboratory; not in scope of accreditation)

Exeter Analytical 440 Elemental Analyzer

Combustion method

%	Theoretical	Result 1	Result 2	Mean
C	42.86	43.00	42.87	42.94
H	3.60	3.47	3.40	3.44
N	24.99	24.88	24.83	24.86

**SIGMA-ALDRICH®**



### HOMOGENEITY ASSESSMENT

Homogeneity was assessed in accordance with ISO Guide 35. Completed units were sampled using a random stratified sampling protocol. The results of chemical analysis were then compared by Single Factor Analysis of Variance (ANOVA). The uncertainty due to homogeneity was derived from the ANOVA. Heterogeneity was not detected under the conditions of the ANOVA.

Analytical Method: HPLC

Sample size: ~20mg

### UNCERTAINTY STATEMENT

Uncertainty values in this document are expressed as Expanded Uncertainty ( $U_{\text{cm}}$ ) corresponding to the 95% confidence interval.  $U_{\text{cm}}$  is derived from the combined standard uncertainty multiplied by the coverage factor  $k$ , which is obtained from a  $t$ -distribution and degrees of freedom. The components of combined standard uncertainty include the uncertainties due to characterization, homogeneity, long term stability, and short term stability (transport). The components due to stability are generally considered to be negligible unless otherwise indicated by stability studies.

### STABILITY ASSESSMENT

Significance of the stability assessment will be demonstrated if the analytical result of the study and the range of values represented by the Expanded Uncertainty do not overlap the result of the original assay and the range of its values represented by the Expanded Uncertainty. The method employed will usually be the same method used to characterize the assay value in the initial evaluation.

Long Term Stability Evaluation - An assessment, or re-test, versus a Compendial Reference Standard may be scheduled, within the 3 year anniversary date of a release of a Secondary Standard. The re-test interval will be determined on a case-by-case basis.

Short Term Stability Study - It is useful to assess stability under reasonably anticipated, short term transport conditions by simulating exposure of the product to humidity and temperature stress. This type of study is conducted under controlled conditions of elevated temperature and humidity.



QC Supervisor



QA Supervisor

### APPENDIX

Original Release Date:	06 October 2016
Stability Test Date:	30 April 2017
Requalification Test Date:	30 April 2017

